BG: ИНСТРУКЦИЯ ЗА УПОТРЕБА
EN: MANUAL INSTRUCTION FOR USE
DE: MANUELLE BEDIENUNGSANLEITUNG
GR: ΟΔΗΓΙΕΣ ΧΕΙΡΙΣΜΟΥ ΧΕΙΡΙΣΜΟΥ
NL: GEBRUIKSAANWIJZING
HU: HASZNÁLATI UTASÍTÁSI UTASÍTÁSOK
RO: INSTRUCȚIUNI DE UTILIZARE

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EN

**MANUFACTURER KEEPS THE RIGHT TO CHANGE SPECIFICATIONS OF THE PRODUCT WITHOUT PRE-WARNING!
IMPORTANT! READ CAREFULLY AND KEEP FOR FUTURE REFERENCE!**

RULES AND RECOMMENDATIONS FOR SAFETY USAGE

- Follow the advises for maintenance, described in the manual.
- Device can be used for professional purposes or for personal needs.
- The device must be used for purposes, described in this manual.
- This device must be used only in temperature range of the environment, from 10 or 40 degrees.
- The device is not intended for using by children. Medical products are not toys!
- The device is constructed for practical application but it can not replace a visit to the doctor.
- Children must not play with the device.
- Cleaning and service can not be performed by children.
- Storage the device away from children. There is swallowing danger! Children can swallow the packing material.
- The device must always storage in clean and dry environment.
- Do not expose the thermometer on extreme temperature conditions, bigger than 55° C and less than- 20° C .

- Do not use the device in relative humidity, higher than 85%.
- Attention! The safety glass over the lens is the most fragile part of the thermometer.
- Do not touch the glass of the infrared lens with fingers.
- Clean the glass with cotton, lightly moisture with 95° degrees alcohol towel.
- Do not expose the thermometer on the harmful impact of sunlight or water.
- Do not allow access for water or other liquids to the thermometer to avoid short circuit.
- Do not place and do not storage this product on place, from where it can fall or to be pushed in sink or bath tub.
- Do not place or drop in water or other liquid.
- Never drop the device.
- Use this product only as intended, as it is described in this manual. Do not make improvements or changes over the product.
- This can lead to its incorrect work or to harming the user.
- Before every usage, check the product for traces of damage. Do not use it if you find damage on any part.
- If a problem with the device appear, contact with the trader, from which you have bought it. Do not try to repair the device by yourself. Otherwise the warranty will be canceled.

➤ **IMPORTANT BATTERY GUIDELINES:**

1. An adult must work with the battery.
2. Do not allow children to play with the batteries. Swallowing hazard! Little children can swallow the batteries and to suffocate. Because of this, storage the batteries on unavailable for children place!
3. Always remove exhausted batteries and do not throw them together with the household waste, but on the intended places. They are recyclable.
4. Use batteries only AAA size.
5. Alkaline batteries are recommended.
6. When you change with new batteries, always change all batteries. Children must not attend by the time of changing the batteries.
7. Always place the batteries correctly according their polarity (+/-), marked on the bottom of batteries. If they are not correctly places, this can lead to short circuit or electric shock.
8. Do not mix old with new batteries.
9. Do not mix alkaline, standard (Carbon-Zinc) or rechargeable (Nickel- Cadmium) batteries.
10. Do not put the batteries in fire, because they can explode or to leak.
11. Keep the batteries from direct heating sources.
12. If you will not use the device for long time, take the batteries out from the battery compartment. This way you will prevent the device from battery leaking or corrosion.
13. If when changing the batteries you establish that the batteries leak, put yourself safety gloves and clean the battery compartment with dry towel.
14. Non-rechargeable batteries must not be recharge.
15. Do not disassemble, open or shatter the batteries.

DEVICE PURPOSE

The device is infrared thermometer, intended to measure the temperature on babies' and old people's forehead, without having any contact to human body. It can be used by the users in home or hospital environment.

INTRODUCTION

Contactless thermometer is made by using the newest infrared technology. This allows measurement of the temperature on the superficial temporal artery (TA) to be made from distance from the forehead from 3 to 5 cm. Precise, fast and wireless, this thermometer is the most suitable for measuring the temperature without risk at the moment. Following studies prove that the method of measuring temperature of TA is more precise than using tympanic thermometry and better tolerated than rectal thermometry (1).

However, as other types of thermometer, Easy Check must be used correctly to obtain reliable and stable results. Therefore we advise carefully to read the manual and the safety measures before using the product.

(1)Greenes D., Fleisher G., Accuracy of Noninvasive Temporal Artery Thermometer for Use in Infants . Arch Pediatr Adolesc Med 2001;155:376.

PRECAUTIONS FOR USE

JXB-183 pre-set at factory. It is not necessary to calibrate the device when starting it up.

In order to obtain reliable and stable results, we recommend you, each time, when there is significant change in the ambient temperature, to allow to Easy Check to acclimatize to this ambient temperature from 15 to 20 minutes before using it. It is important to leave interval from 3-5 second between measurements.

All objects, solid, liquid or gas emit energy by radiation. The intensity of this energy depends of the temperature in the object. The infrared thermometer Easy Check is therefore able to measure the human's temperature, only by the energy that person emits. This measurement is possible, thanks to anexternal temperature probe on the device, which permanently analyses and register the temperature of the environment.

Therefore, immediately after the thermometer is closed to the body and the radiation sensor is activated, the temperature is immediately measured through detecting infrared heat, caused by arterial blood flow. The temperature of the body can be measured also without any interference from the heating of the surrounding environment.

DIFFERENT METHODS FOR TEMPERATURE MEASUREMENT

A. Internal temperature

The measurement of the internal temperature is the biggest method for body temperature readings and this includes measurement of the pulmonary artery, using a catheter equipped with a thermal probe that can measure the temperature on site. The same method is used for probes, measuring the temperatures in the esophagus. Although, such an invasive method requires specific equipment and expertise.

B. Rectal thermometry

The rectal temperature is changing slowly, compare to the internal body temperature. It is proved that the rectal temperature remains increased for long time after the internal temperature of the patient begin to drop and in vice versa. Furthermore, rectal perforations occur in a result of this method and without sterilization techniques, rectal thermometer can spread germs, often found in faeces.

C. Oral thermometry

Oral temperature is under influence of the recent ingestion of food and drinks and from mouth breathing. For measuring the oral temperature, the mouth must be closed and the tongue to be relaxed for 3-4 minutes, which is hard difficult task for children.

D. Axillary (armpit) temperature

Although that the armpit temperature can be easily measured, it is proved that it does not show exact measurement of the children internal temperature. For measuring this way, the thermometer must be tightened hard over the armpit artery. Although lower sensitivity and comparative inaccuracy of the temperature of the armpit for detection of fever, this method is recommended by the American Academy of Pediatrics as screening test for fever in newborns.

E. Tympanic thermometry

In order to obtain precise temperature reading, good command of the measurement technique is required. The thermometer probe must be placed as close as the warmest part of the external ear canal.

NORMAL TEMPERATURES ACCORDING TO THE MEASURING METHOD	
MEASURING METHOD	NORMAL TEMPERATURE
RECTAL	36.6 °C – 38 °C
ORAL	35.5 °C – 37.5 °C
AXILLARY	34.7 °C – 37.3 °C
AURICULAY	35.8 °C – 38 °C
TEMPORAL	35.8 °C – 37.8 °C

Human body temperature varies through the whole day. Also it can be influenced by many external factors: age, sex, type and thickness of the skin.

➤ **Advantages of temporal artery (TA) temperature**

Infrared arterial temperature can be measured, using device, placed on the forehead in the area of the temporal artery. It is proved that this relatively new method for measuring the temperature is more precise than using of ear thermometers and better tolerated than the rectal method.

The thermometer Easy Check, item number JXB-813, is intended to show immediately the forehead temperature without any contact with the temporal artery. Because temporal; artery is too close to the surface of the skin and therefore accessible and considering that the blood flow is permanent and regular, this allows precise temperature measurement. This artery is connected with the heart through the carotid artery, which is directly connected with the aorta. This form part of the main part of the main stem of the arterial system. The effectivity, speed and comfort from measurement of the temperature from this part of the body, make this method ideal in comparison with other methods for temperature measurement.

NORMAL TEMPERATURE ACCORDING AGE

Age	°C	°F
0-2 years old.	36.4 – 38.0	97.5 – 100.4
3-10 years old.	36.1 – 37.8	97.0 – 100.0
11-65years old.	35.9 – 37.6	96.6 – 99.7
> 65 years old.	35.8 – 37.5	96.4 – 99.5

Body temperature of infants and small children is higher and decrease with the age. the children quickly and often experience more severe temperature fluctuations, caused by leaps in growth.

➤ **Practical considerations when taking a temperature**

1. To be sure that the measurement is precise and accurate, it is important every user to receive adequate information and training for temperature measurement with such device.
2. It is important to remember that, through procedures as temperature measurement to seem easy, they should not be underestimated.
3. The temperature must be measured in neutral context. The patient must not be done any energetic physical activity, before his temperature to be measured and the room temperature must be moderate.
4. Be careful from physiological variations in temperature, which must be taken under attention in the measurement results: temperature increase with 0,5 degrees between 6 a.m. and 3 p.m. Women have about 0,2 degrees higher temperature. Their temperature varies according their period. In the second half of their period, the temperature increase with 0,5 degrees, as in early stage of pregnancy.
5. When sitting, the temperature is lower with about 0,3 to 0,4 degrees, than when we stand upright.
6. The temperature is affected from different factors, as individual exchange of substances; age; clothing; environment temperature; previous body activity and to a lesser extent from mental activity; place, on which it is measured.

➤ **How to measure temperature**

Hold the thermometer against forehead, in the right area, on distance from 3 to 5 cm. Press the button for temperature measurement and it will immediately appear on the display.

ATTENTION! Reliability of measurement can not be guaranteed if the temperature is measured on other part of the body (for example hand, torso...).

➤ **Restrictions**

Please, keep the following restrictions before measuring the temperature, for ensuring real and reliable result:

- Remove any hair from the forehead.
- Wipe the forehead, if sweating appeared.
- Avoid air flows (for example. air conditioners and other)
- Leave from 3 to 5 second interval between 2 measurements.
- With every significant change in the room temperature, leave Easy Check to adapt to this temperature at least 15 minutes before using.

DEVICE DESCRIPTION

Look page 2

Attached part type BF: Sensor

1. Infrared lenses; 2. LCD display; 3. Alarm; 4. Memory; 5. Celsius or Fahrenheit; 6. Mode; 7. On/Scan button

CHARACTERISTICS

- 1.Special design for measuring body temperature from distance 3-5 cm from forehead.
- 2.Reliable and true measurement, thanks to the advantage “Infrared reporting system”.
- 3.Sound alarm, if the temperature is more than 38 °C .
- 4.Remembers last 32 temperature measurements.
- 5.Blacklits LCD digital display screen.
- 6.The temperature is recorded as in Celsius, both in Fahrenheit .
- 7.Automatically turning off (≤30 secs.), for energy saving.
- 8.Long time for using (100 000 measurements).
- 9.Practical and easy for using.

➤ Additional usage:

Easy Check can be used both for temperature measuring in baby feeding bottle and temperature of bath tub water (using mode for surface temperature measuring) or room temperature (using room mode).

INSTRUCTIONS FOR USING

1. Put the batteries in the compartment correctly according its polarity.
2. While using this for the first time or when putting new battery, pls wait for 10-15 minutes for heating the device. This will allow to the device to acclimate to the room temperature.
3. Press the on/scan button, point the device to the forehead (look diagram for pointing the thermometer), on distance from 3 to 5 cm. press the button ON/ Scan standby mode and the measuring will finish then the temperature shows itself on the display or a sound is issued.
4. Measuring time is one second. **Note: Do not change the position of the thermometer before measuring to be finished.**
5. Before temperature measuring, ensure that you have removed the hair or sweating from the forehead.

MENU SETTINGS AND FUNCTIONS

1. Turning the device on:

Press the ON/ Scan button, after second, the panel in the display will be in standby mode with signs „----°C“ и „---°F“. Then press the ON/SCAN button once more and the result will show up automatically after 1 second. If you do not use the thermometer anymore, it will turn off automatically after 30 seconds.

2. Mode switching

- a. Press the button “MODE” and the display will show: Body ...°C
- b. Press the button “MODE” again and display will show : Room ...°C
- c. Press again the button “MODE” and display will show: Surface Temp....°C

Note: The thermometer is set by default on „BODY“ mode.

IMPORTANT!

Surface temperature differs from the internal temperature of human body. For measuring the body temperature, always use „BODY“ mode. Please, ensure that you have choose “SURFACE TEMP” mode only for reporting the external area.

3. In turned on state, press the button, °C/°F“, to switch between Celsius and Fahrenheit.
4. In turned on state, press the button “MEM” (Memory), which will show the last measured temperature and will show last 32 measurements.
If again, in turned on state, you press the button “MEM” and hold it for 5 seconds, all data will be deleted. After this press again „MEM“ button and the display will show „CLR“.
5. In tuned on state, press button „“, to activate or deactivate the alarm. When the display shows “ON”, this means the alarm is activated. When the display shows “OFF”, this means that the alarm is deactivated.
6. Device calibration though menu F4

When there is difference between Easy Check and mercury thermometer, you believe that the mercury thermometer is more accurate but not comfortable for using. You can use calibration function to tune Easy Check and to receive the value, same as this on the mercury thermometer. Furthermore, when you use Easy Check for people with different skin color (for example from yellow, white or black race), you can also use calibration.


➤ Calibration instructions

1. In turned on state, press the buttons “MODE” and „°C/°F“ at the same time for 3 seconds. The display will show: “F4”.
2. Press button “MODE”, to increase with 0.1 °C, press „°C/°F“, to decrease with 0.1 °C.
3. Press “MEM” button, to keep the settings.
4. In case of seasonal or ecological changes, a check and device regulation.

➤ Battery change:

Display: when LCD display shows this symbol „“, it means batteries are exhausted.

Open the cap and replace batteries, as you be careful to keep the right polarity, when placing new batteries. An error with this action can lead to damage of the device and warranty canceling. Never use rechargeable batteries. Use only single use batteries.

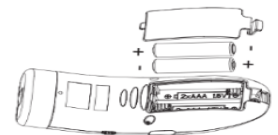TECHNICAL CHARACTERISTICS

1. Normal using conditions
Environment temperature: 10 °C – 40 °C (50 °F – 104 °F)
Relative humidity: ≤ 85%
2. Batteries: DC 3V (2 бр. AAA batteries)
3. Product size: 130 x 45 x 55 мм (L x W x H)

4. Product weight (with batteries): 67 gr.
5. Display temperature resolution: 0.1°C (0.1°F)
6. Measuring range:
 „BODY“ mode : 32°C – 42.9°C (89.6°F – 109.2°F)
 „Surface temperature“ mode : 0°C – 60°C (32 °F - 140 °F)
 „Room“ mode : 0°C - 40°C (32 °F - 104 °F)

7. Precision

- 32.0 °C - 34.9 °C (89.6 °F – 94.8 °F) ± 0.3 °C (± 0.6 °F)
 35.0 °C - 42.0 °C (95 °F – 107.6 °F) ± 0.2 °C (± 0.4 °F)
 42.1 °C - 42.9 °C (107.8 °F – 109.2 °F) ± 0.3 °C (± 0.6 °F)

8. Power: ≤300mW

9. Accuracy: ± 0.3 °C (0.6 °F)

10. Measuring distance: 3cm – 5cm (1.2in – 2in)

11. Automatically turning off: ≤30 secs

12. Memory: 32 measurements

• Note: Contactless infrared thermometer, model Easy Check, item number JXB-183 can take temperature readings under 32°C and over 42,9°C (89,6°F до 109,2°F), but the precision is not guaranteed out of this range.

➤ Durability of the product

Easy Check is made for intensive and professional usage, its durability is guaranteed for 100 000 measurements.

PRODUCT MAINTENANCE

1. Safety glass over the lens is easily broken part from the thermometer, please, be careful with it.
2. Clean the glass with cotton towel, lightly moisture with 95° alcohol.
3. Do not expose the thermometer on sunlight or water. Their influence will damage the product.
4. Clean the whole device with soft cloth. Do not wash it under running water! Do not use aggressive cleaning detergents, as bleach, those with abrasive particles or on alcohol basis.
5. If you will not use the device for long time, take the batteries out, to avoid leaking or corrosion.
6. Do not storage the device on high or low temperatures or high humidity (look technical characteristics), on direct sunlight, close to direct heating sources or dusty places. Otherwise measurement deviations can be caused.

ACCESSORIES

- Manual 1 pcs.
- AAA alkaline batteries 2 pcs.
- Storage case 1 pcs.

GUIDELINES

This device corresponds to the European directive 93/42/EEC, related to medical products, ISO 80601-2-56 and European standard EN60601-1-2 and it is object of concrete safety equipment, regarding electromagnetic compatibility.

TROUBLESHOOTING

If you have problems, while using the thermometer, please, take a look at the guidelines in this manual for solving the problem. If the problem still exists, contact with the trade agent, which you bought the product from.

1. Display shows temperature, higher than 42,9 °C (109,2 °F)

The temperature is in Fahrenheit. Change the measurement to Celsius.

2. Display shows temperature lower than 32 °C (89,6 °F)



For measuring the surface temperature, press button „BODY“ and set for measurement on „BODY“. If the device is in mode „Surface temperature“, the temperature 32 °C (89,6 °F), which display will show, taking into account the external body temperature, instead internal.

3. Display show the message “Hi”

When using contactless thermometer, model Easy Check, the message “Hi” can show on the display. In this case the temperature is either over the chosen range or over 42,9 °C (109,2 °F) in „BODY“ mode.

4. Display shows message “Lo”

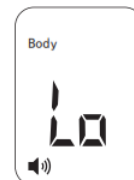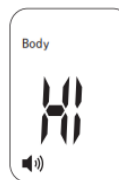When using the contactless thermometer, model Easy Check, the message “Lo” can show on the display. In this case, recorded temperature is either under chosen range or under 32 °C (89,6 °F) in „BODY“ mode.

The display shows this message because of many reasons. Please, take a look at the above list with the main possible problems:

Reasons, the display shows the message “Lo”	Advise
Temperature measurement is difficult from hair or sweating.	Ensure that there are no obstacles or moisture before temperature measurement.
Temperature measurement is difficult from air flows or sharp changes in environment temperature.	Ensure that there are no air flows around the area, in which the temperature is measured, this can affect the infrared reading.
Successive measurements of the temperature are too close one after another and the thermometer do not have time to restart.	Leave interval 3-5 secs. Between measurements; recommended is pause for 15 sec.
Measurement area is away from the thermometer.	Measure the temperature on recommended distance (around 3-5 cm, 1.2 ~2.0 inches).

SYMBOL EXPLANATION

Symbol	Explanation
	Trade mark.
	IEC 60417-5333, type attached part BF
	IEC 60417-5031, Direct current
IP22	Protected access with fingers to dangerous parts and against vertical falling of water drops, when it is tilted up to 15 °C
	Directing to the instruction/Manual.
	DISPOSING: Do not throw this product with household waste. Collect separately, because special disposal of those kind of waste is required.
SN	Specific serial number.

DECLARATION OF CONFORMITY

Manufacturer guide and declaration – electronic immunity

Easy Check is intended for using in electromagnetic environment, specified below. The client or user of JXB-183 must be sure that the product is used in such environment.

Immunity test	IEC 60601 level test	Level of compliance	Electromagnetic environment - guide
Electromagnetic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors must be wooden, concrete or ceramic tiles. If the floors are covered with synthetic materials, relative humidity must be at least 30 %.
Electrical fast transient, explosive IEC 61000-4-4	± 2 kV 3a power lines ± 1 kV for input-output lines	NOT APPLICABLE	The quality of the basis power must be this of typical trade or hospital environment.

Manufacturer guide and declaration – electronic immunity

Easy Check is intended for using in electromagnetic environment, specified below. The client or user of JXB-183 must be sure that the product is used in such environment.

Immunity test	IEC 60601 level test	Level of compliance	Electromagnetic environment - guide
Waves IEC 61000-4-5	± 1 kV line to line ± 2 kV line to ground	NOT APPLICABLE	Quantity of the basis power must be this on typical trade and hospital environment.
Voltage dips, short interruptions and voltage variation of the input power lines IEC 61000-4-11	< 5% UT (>95% dip in UT) for 0,5 cycle 40% UT (60% dip in UT) for 5 cycle 70% UT (30 % dip in UT) for 25 cycle <5% UT (>95% dip in UT) for 5 secs	NOT APPLICABLE	Mains power quality should be that of a typical commercial or hospital environment. If the user of the “JXB-183” requires continued operation during power mains interruptions, it is recommended that the “JXB-183” be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) Magnetic field IEC 61000-4-8	3 A/m	3 A/m	Frequency magnetic fields must be on the levels of the typical location of the typical trade or hospital environment.

NOTE: UT is pressure of the alternating voltage before applying the testing level.

Guides and declarations from the manufacturer – electric immunity

Easy Check is intended for using in electromagnetic environment, aspecified below. The client or user of JXB-183 must ensure, that it is used in such environment.

Immunity test	IEC 60601 level test	Compliance level	Electromagnetic environment- guidelines
			Portable and mobile radiofrequency equipment must not be used close to none of the parts of “Easy Check”, including cable, except the recommended distance, calculated from the equation, applicable for transmitter’s frequency.
Conducted radiofrequency IEC 61000-4-6	3 Vrms 150 kHz for 80 MHz	NOT APPLICABLE	$d = \left[\frac{3E}{V} \right] \sqrt{P}$ $d = \left[\frac{3E}{E} \right] \sqrt{P}$ $d = \left[\frac{7}{E} \right] \sqrt{P}$
			80MHz for 800 MHz 800MHz for 2.5 GHz

Immunity test	IEC 60601 level test	Compliance level	Electromagnetic environment- guidelines
Broadcasted radio frequency IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	Where P is the maximum nominal power in wats (W), according to the manufacturer’s transmitter and “d” is the recommended distance in meters (m). Power of the filed from fixed RF transmitters, as described from survey of electromagnetic objects, “a” must be less than level of compliance in every frequency range “b”.

Immunity test	IEC 60601 level test	Compliance level	Electromagnetic environment- guidelines
NOTE 1: For 80 MHz and 800 MHz is applicable the highest frequency range. NOTE 2: These guidelines may not apply in all situations. Electromagnetic multiplication is not influenced by absorption and reflection from structure, objects and people.			
a	<p>The strength of the field from stationary transmitters, as main situations for radio (cellular/wireless) telephones and terrestrial mobile radios, amateur radio, AM and FM broadcasting and TV broadcasting, theoretically it can not be provided with accuracy. To evaluate the electromagnetic environment, because of fixed radiofrequency transmitters, it have to be provided survey of electromagnetic objects. If the measured power of the field, on which JXB-183 is used, increases the applicable compliance level with RF above, medical JXB-183 must be monitored, to check the normal working. If observed abnormal performance, additional measures must be taken, as reorientation or moving of JXB-183.</p> 		
b	At frequency, higher than the range 150 KHz to 80 MHz, the strength of the field must be less than 3 V/m.		
Recommended distances for separating between portable and mobile RF equipment and the medical "JXB-183"			
"JXB-183" is intended for using in electromagnetic field, in which the emitted RF interruptions are controlled. The client or user of the medical "JXB-183" can help for preventing the electromagnetic interruptions, as it keeps to the minimal distances between portable and mobile RF communication equipment (transmitters) and "JXB-183", as the recommended below, according maximal power of the communication equipment.			
The calculated maximum output power of the transmitter	Separation distance, according the frequency of transmitter.		
	150 KHz to 80 MHz $d = [3,5/V1] VP$	80 MHz to 800 MHz $d = [3,5/E1] VP$	800 MHz to 2,5 GHz $d = [7/E1] VP$
0,01	/	0,12	0,23
0,1	/	0,38	0,73
1	/	1,2	2,3
10	/	3,8	7,3
For transmitters, calculated with maximal output power, which are not recorded above, the recommended separating distances d in meters (m) can be predicted, using the equal, applicable to the transmitter's frequency, where P is the maximum output power, measured in wats (W), according to the manufacturer's transmitter. NOTE 1: For 80 MHz and 800 MHz, the separating distance for the higher frequency range applies. NOTE 2: This guideline may not apply in all situations. The electromagnetic multiplication is influenced by absorption from structures, objects and people.			

Produced for CANGAROO in China

Manufacturer: Guangzhou Jinxinbao Electronic Co., Ltd.

Address: 38 Huanzhen Xi Road Dagang Town, Nansha, Guangzhou, China

Importer: Moni Trade Ltd.

Address: Bulgaria, Sofia, Trebich quarter, ul. Dolo 1

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Producător: Guangzhou Jinxinbao Electronic Co., Ltd.
IMPORTATOR : MONI TRADE LTD