Tender Care

Digital Soother Thermometer

MANUAL INSTRUCTION

designed in EU

www.lorelli.eu

V1.1W

find us on

CONTENT

EN	Manual Instruction.....	3
BG	Инструкция за употреба.....	4
GR	ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ ΚΑΡΕΚΛΑΣ.....	5
ES	Instrucciones de uso.....	6
RU	Инструкция по эксплуатации.....	7
HU	Bedienungsanleitung.....	8
RO	Instrucțiuni de utilizare.....	9
IT	Istruzione per l'uso.....	10
FR	Mode d'emploi.....	11
AR	دليل دليلاتعليمات.....	12
DE	Bedienungsanleitung.....	13
CZ	Návod k použití.....	14
SRB HR ME BIH	Uputstvo za upotrebu.....	15
NL	Gebruikshandleiding.....	16
AL	Insruksion për përdorim.....	17
TR	Kullanım talimatı.....	18
PL	Instrukcja użytkownika.....	19
MK	Упатствата за употреба.....	20

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EN OPERATING INSTRUCTIONS

Dear Customer, Thank you for purchasing Digital Soother Thermometer. Before using this product, please read this instruction carefully.

1. Description

2. Feature

- LCD screen display.
- Memory function.
- Automatically taking the temperature and shut down automatically.
- No mercurial pollution or safe concern.

3. Specification

- Measurement Scope: 32.0 °C - 42.9 °C (89.6 °F - 109.2 °F)
- Measurement Time: Three minutes +/- One minute
- Temperature Accuracy: +/- 0.1 °C (32.2 °F)
- Battery: 1.50 VDG AG3
- Display Mode: LCD (liquid crystal display)
- Memory: Display the previous temperature.
- Timer: When the measure completely, it will ring 10 times.
- Dimension: 55 x 50 x 50 mm
- Weight: About 16g

4. How to use it?

■ Press on/off button, wait for LCD display 188.8, to make sure the thermometer working well and battery stay good condition, then it will display the previous saving memory as well as activate the working figure 37.0 °C (LO °F) Flash °C (°F), if the environmental temperature is higher than 32.0 °C (89.6 °F) it will directly display the figure and begin measuring the temperature.

- If the temperature is higher than 42.9 °C (109.2 °F) the LCD screen will display H °C (H °F). If the temperature is below than 32.0 °C (89.6 °F) screen will display LO °C (LO °F).
- When the temperature is complete LCD screen will display °C (°F) and stop flashing with 10 times ring music.
- When you have got the temperature result, please make sure shut down the thermometer in order to save battery's life. If not, operation system will automatically turn off the thermometer in 10 minutes.
- If the measurement is shut down, please repeat the above procedure to begin the temperature measurement again.

5. Precautions

- Before measuring the temperature, please check the soother first. If it has broken problem, please do not use it.
- Do not put the thermometer in high heat or expose to direct sunshine and prevent touching any chemical solvent to avoid any unexpected situation happen.
- Please, do not try to disassemble, repair or re-build the thermometer.
- Digital Soother Thermometer is an electronic measurement. Please, keep away from your child.
- Digital Soother Thermometer is for mouth measuring only.
- The battery is sealed installed. Please, do not take part it from the thermometer by yourself.

6. Clean the thermometer

- Put the soother into the alcohol solvent in one second. LCD display "Screen can only be cleaned by fabric or by clean water".
- This thermometer has waterproof function. You can put the thermometer in clean water for simply washing.

Caution: This thermometer can not be put into boiled water.

7. Limited warranty

The warranty does not cover the damage of improper use or the battery running out. If the unit does not function properly due to defective parts or assembly we will repair it free of charge or replace with a new one.

Symbol	Description
	Type BF applied part.
	Information about a manufacturer
	Please read the instructions carefully.
IP22	Rate degree of protection against access to water and inhomogeneous substances
	Waste electrical materials should be sent to a dedicated collection point for recycling.
	EU safety requirements
	European authorized representative

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Authorized Signature(s)

Declaration of Compliance

Company Name, Invoice Company Address: FUDAKANG INDUSTRIAL CO. LIMITED,
Room 2702,27/F., Omega Plaza,32 Dundas Street, Mongkok, KLN Hong Kong, China

Company Name, Factory Address: FUDAKANG INDUSTRIAL CO. LIMITED, 301, 3rd
Floor, Block 1, No.60 Shangshuxia Road, Changping Town, Dongguan, Guangdong,
P.R. China

Product Item No: BT-A51

Description : DIGITAL SOOTHER THERMOMETER

For our model BT-A51, delivered to DIDIS Ltd as 1025006 DIGITAL SOOTHER
THERMOMETER:

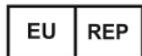
1. We hereby declare that the above model is compliant to the Restriction of Hazardous Substances in Electrical and Electronic Equipment Directive (RoHS) 2011/65/EU.
2. We hereby agree that only RoHS compliant materials and components will be used in the assembly of the above model. Under no circumstances will non-compliant materials and components be used in place of approved compliant equivalents during the life of the product.
3. This device corresponds to the European directive 93/42/EEC relating to medical devices, ISO 80601-2-56 and European standard EN60601-1-2 and it is subject to specific precautions, regarding the electromagnetic compatibility.

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**MedPath GmbH, Add: Mies-van-der-Rohe-Strasse 8,
Importator:Didis Ltd.